

Sheep Ear Tagging

1. Preparation of Dual Tags
 - a. Separate the EID and VID tags.



- b. Remove the connectors by twisting.



** For cold weather, it is advisable to keep the tags in a warm container.


2. Insertion of Tags to the Applicator
 - a. Place the Female Part (Socket) of the Tag to the Holed Jaw.



- b. Place the Male Part (Pin) of the Tag to the Narrow Jaw.



- c. Slightly grip the Applicator to check if the Pin will go to the Socket.

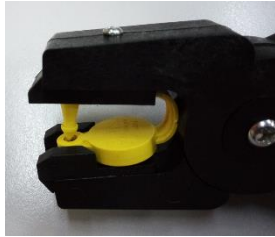



3. Apply disinfectant onto the Tag and Applicator.

4. Tagging the Sheep


- a. Firmly hold the head and the ear of the Sheep.
b. Position the Female part (Socket of the Tag inside the ear.



- c. The tag must be positioned 1/3 from the base of the ear and must not hit the cartilage or blood vessels (vein) of the ear.

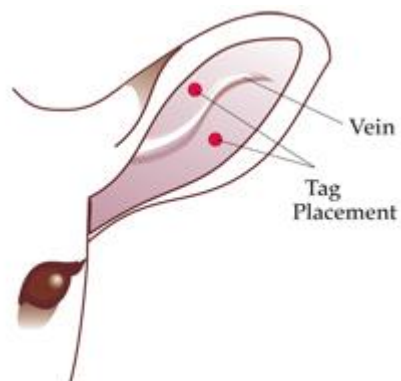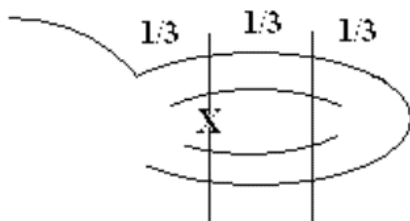

d. For lambs, allocate at least 5 mm gap from the edge of the ear to allow growth.

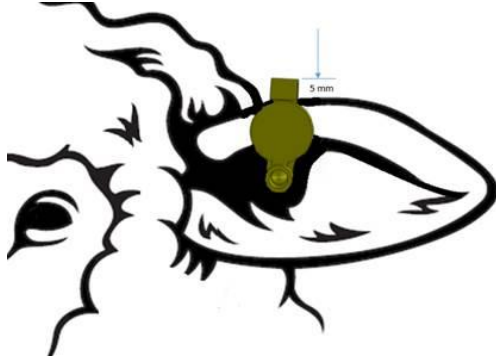

d. For older sheep or sheep with thick ears, it is advisable to apply the tags to the back of the ear.



- e. Position the Male Part (Pin) of the Tag to be almost **perpendicular** to the ear.
- f. Squeeze the applicator quickly and firmly until the Pin goes through Ear and the Socket.

To view the video of how to install the Tag, click on the below link :

For further enquiries or questions, please kindly email us at : enquiries@ptstechnologies.com and we will channel you to our renown Distributors near you (Distributors: Ireland, UK, Romania, Germany, New Zealand, Australia, USA, Thailand, Indonesia, Vietnam).